

BioClean-D™ Coverall with Hood and Integrated Boots BDFC

Hooded non-sterile coverall with boots, providing chemical and liquid splash protection

- Reduced contamination risks: The BioClean-D Coverall with Hood and Integrated Boots BDFC is low-linting and ISO cleanroom-compliant thanks to its lightweight CleanTough™ composition, ensuring contamination risks are kept to a minimum
- Optimized protection: This non-sterile disposable coverall features a front zip with a protective flap, which guards against liquid splashes and chemicals, while thumbs loops and an elasticated hood, back, cuffs and ankles offer a secure fit
- Enhanced features: This antistatic coverall is also equipped with integrated boots and slip-resistant soles, for added workplace safety in laboratory environments

Key Features and Benefits

- CleanTough™ fabric: For lightweight, low-linting body protection
- Front zip and protective flap:
 Contaminant and liquid splash barrier
- Integrated boots and anti-slip soles: For head-to-toe cleanroom PPE

Industries

- Controlled and Critical Environments
- Production and Manufacturing
- Pharmaceutical Manufacturing
- Biotechnology Manufacturing
- Medical Device Manufacturing









BioClean-D™ Coverall with Hood and Integrated Boots BDFC

TECHNICAL DATA SHEET

PRODUCT INFORMATION

Material	CleanTough™
Audit Standards	ISO 14001, Manufacturing QMS Audit Standards ISO 9001, PPE Regulation 2016 425 Module D, NEBB Certified Cleanrooms
Standards	ASTM F739, CE 0598, EN ISO 13688:2013, EN 1149-5:2018, EN 13934-1, EN 13935-2, EN 530, EN 6530, EN 7854, EN 863, EN 9073-4, EN ISO 14325, Category III, EN ISO 13982-1:2004 + A1:2010, EN 13034:2005 + A1:2009
Packaging Overview	One piece aseptically folded per sealed inner PE bag; one inner bag per sealed outer PE bag; 20 outer bags per lined carton (20 pieces). Please note: sizes XXXL & XXXXL 15 pieces per carton. All sizes subject to minimum order quantity and lead time.
Country Of Origin	China
Cleanroom Class	Class 10/ISO 4
Shelf Life	Five (5) years from date of manufacture.
Construction	Bound seams with single needle stitching
Characteristics	*NOTE: BioClean CleanTough material is static dissipative and, with a charge half decay time of 0.07 sec, and so are ideal for use in a static-safe environment.

PARTICLE SHEDDING TEST RESULTS

TEST	RESULT	
Particle Shedding (Helmke Drum Test)	≥ 0.5Qm (counts/min) <2000	

ASTM F739-12 TEST METHOD RESULTS

DRUG	Mean Breakthrough Time (MBT), Minutes Breakthrough of the test chemical is deemed to have occurred when the permeation rate has reached 0.1 Qg/cm² /min
CISPLATIN	>240
CARMUSTINE	<6
CYCLOPHOSHAMIDE	217 (275,162,215)
DOXORUBICINHYDROCHLORIDE	>240
5-FLUOROURACIL	>240
METHOTREXATE	>240
ETOPOSIDE	>240
PACLITAXEL	<10
THIOTEPA	30 (28,30,33)

Results achieved under controlled laboratory conditions, by accredited external testing laboratory.

SIZE CHART

BDFC-S; Size: S, Chest: 84-92cm (33"-36"), Height: 164-170cm (5'4"-5'6") BDFC-M; Size: M, Chest: 92-100cm (36"-39"), Height 170-176cm (5'6"-5'9") BDFC-L; Size: L, Chest: 100-108cm (39"-42"), Height: 176-182cm (5'9"-6'0") BDFC-XL; Size: XL, Chest: 108-116cm (42"-45"), Height: 182-188cm (6'0"-6'2") BDFC-2XL; Size: 2XL, Chest: 116-124cm (45"-48"), Height: 188-194cm (6'2"-6'4") BDFC-3 XL; Size: 3XL, Chest: 124-132cm (48"-51"), Height: 194-200cm (6'4"-6'6") BDFC-4XL; Size: 4XL, Chest: 132-140cm (51"-54"), Height: 200-206cm (6'6"-6'8")





BioClean-D™ Coverall with Hood and Integrated Boots BDFC

MATERIAL PERFORMANCE TEST RESULTS

TEST	RESULT	PERFORMANCE CLASS
Abrasion Resistance	10 to 100 cycles	1
Flex Cracking Resistance	5,000 to 15,000 cycles	3
Puncture Resistance	8 N	1
Trapezoidal Tear Resistance Cross Direction (CD)	CD 29.3 N	2
Trapezoidal Tear Resistance Machine Direction (MD)	MD 55.5 N	3
Tensile Strength Cross Direction (CD)	CD 48 N	1
Tensile Strength Machine Direction (MD)	MD 97 N	2
Repellence to Liquids - 30% H ₂ SO ₄	96.3%	3
Repellence to Liquids - 10% NaOH	97.6%	3
Repellence to Liquids - O-Xylene	95.7%	2
Repellence to Liquids - Butan-1-ol	96.6%	3
Penetration by Liquids - 30% H ₂ SO ₄	0%	3
Penetration by Liquids - 10% NaOH	0%	3
Penetration by Liquids - O-Xylene	0%	3
Penetration by Liquids - Butan-1-ol	0%	3
Inward leakage of aerosols of fine particles ¹	Pass	N/A
Resistance to penetration by spray (Spray Test) ¹	Pass	N/A
Seam Strength ²	70 N	2
Electrostatic Charge Half Decay Time, t ₅₀ (secs)	0.07	PASS

^{1.} All interfaces with additional PPE were taped during testing, as required by the standards.

ORDERING INFORMATION

	SIZE	S, M, L, XL, 2XL	3XL, 4XL
BDFC	REORDER NO.	BDFC-S, BDFC-M, BDFC-L, BDFC-XL, BDFC-2XL, BDFC- 3XL	BDFC-3XL, BDFC-4XL

For additional information visit us at www.ansell.com, or call us at

Europe, Middle East & Africa Region

Ansell Healthcare Europe NV T: +32 (0) 2 528 74 00 F: +32 (0) 2 528 74 01

Asia Pacific Region

Ansell Global Trading Center T: +603 8310 6688 F: +603 8310 6699

North America Region

Ansell Healthcare Products LLC US T: +1 800 800 0444 US F: +1 800 800 0445 CA T: +1-800-363-8340

Latin America & Caribbean Region

Ansell Commercial Mexico S.A. de C.V. T: +52 442 248 1544 / 248 3133

Australia

Ansell Limited T: +61 1800 337 041 F: +61 1800 803 578

UK

Ansell Nitritex T: +44 1638 663338 F: +44 1638 668890 **Performance Standards and Regulatory Compliance**












Ansell, $^{\textcircled{8}}$ and $^{\textcircled{m}}$ are trademarks owned by Ansell Limited or one of its affiliates. US Patented and US and non-US Patents Pending: www.ansell.com/patentmarking $^{\textcircled{6}}$ 2024 Ansell Limited. All Rights Reserved.

Neither this document nor any other statement made herein by or on behalf of Ansell should be construed as a warranty of merchantability or that any Ansell product is fit for a particular purpose. Ansell assumes no responsibility for the suitability or adequacy of an end user's selection of gloves for a specific application.

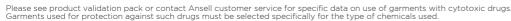



^{2.} Seam not destroyed